

✓ BPL_IEC_validation script with PyFMI

The key library PyFMI is installed.

After the installation a small application BPL_IEC_validation is loaded and run. You can continue with this example if you like.

```
!lsb_release -a # Actual VM Ubuntu version used by Google
```

```
➦ No LSB modules are available.
Distributor ID: Ubuntu
Description:   Ubuntu 22.04.4 LTS
Release:      22.04
Codename:     jammy
```

```
%env PYTHONPATH=
```

```
➦ env: PYTHONPATH=
```

```
!python --version
```

```
➦ Python 3.11.11
```

```
!wget https://repo.anaconda.com/miniconda/Miniconda3-py311_24.11.1-0-Linux-x86_64.sh
```

```
!chmod +x Miniconda3-py311_24.11.1-0-Linux-x86_64.sh
```

```
!bash ./Miniconda3-py311_24.11.1-0-Linux-x86_64.sh -b -f -p /usr/local
```

```
import sys
```

```
sys.path.append('/usr/local/lib/python3.11/site-packages/')
```

```
➦ --2025-02-07 07:03:20-- https://repo.anaconda.com/miniconda/Miniconda3-py311_24.11.1-0-Linux-x86_64.sh
Resolving repo.anaconda.com (repo.anaconda.com)... 104.16.191.158, 104.16.32.241, 2606:4700::6810:20f1, ...
Connecting to repo.anaconda.com (repo.anaconda.com)|104.16.191.158|:443... connected.
HTTP request sent, awaiting response... 200 OK
Length: 145900576 (139M) [application/octet-stream]
Saving to: 'Miniconda3-py311_24.11.1-0-Linux-x86_64.sh'
```

```
Miniconda3-py311_24 100%[=====] 139.14M 196MB/s in 0.7s
```

```
2025-02-07 07:03:20 (196 MB/s) - 'Miniconda3-py311_24.11.1-0-Linux-x86_64.sh' saved [145900576/145900576]
```

```
PREFIX=/usr/local
Unpacking payload ...
```

```
Installing base environment...
```

```
Preparing transaction: ...working... done
Executing transaction: ...working... done
installation finished.
```

```
!conda update -n base -c defaults conda --yes
```

```
➦ Channels:
- defaults
Platform: linux-64
Collecting package metadata (repodata.json): done
Solving environment: done
```

```
## Package Plan ##
```

```
environment location: /usr/local
```

```
added / updated specs:
- conda
```

The following packages will be downloaded:

package	build	
ca-certificates-2024.12.31	h06a4308_0	128 KB
certifi-2025.1.31	py311h06a4308_0	163 KB
Total:		291 KB

The following packages will be UPDATED:

```
ca-certificates      2024.11.26-h06a4308_0 --> 2024.12.31-h06a4308_0
certifi              2024.8.30-py311h06a4308_0 --> 2025.1.31-py311h06a4308_0
```

Downloading and Extracting Packages:

certifi-2025.1.31	163 KB	: 0% 0/1 [00:00<?, ?it/s]
ca-certificates-2024	128 KB	: 0% 0/1 [00:00<?, ?it/s]
certifi-2025.1.31	163 KB	: 100% 1.0/1 [00:00<00:00, 15.59it/s]

Preparing transaction: done
Verifying transaction: done
Executing transaction: done

```
!conda --version  
!python --version
```

🔄 conda 24.11.1
Python 3.11.11

```
!conda config --set channel_priority strict
```

```
!conda install -c conda-forge pyfmi --yes # Install the key package
```

🔄

Preparing transaction: done
Verifying transaction: done
Executing transaction: done

```
# Notes of BPL_IEC_validation
```

Now specific installation and the run simulations. Start with connecting to Github. Then upload the two files:

- FMU - BPL_IEC_Column_system_linux_om_me
- Setup-file - BPL_IEC_explore

```
%%bash
git clone https://github.com/janpeter19/BPL_IEC_validation

Cloning into 'BPL_IEC_validation'...

%cd BPL_IEC_validation

/content/BPL_IEC_validation
```

✓ BPL IEC validation

Author: Jan Peter Axelsson

```
run -i BPL_IEC_explore.py
```

Linux - run FMU pre-compiled OpenModelica

Model for bioreactor has been setup. Key commands:

- par() - change of parameters and initial values
- init() - change initial values only
- simu() - simulate and plot
- newplot() - make a new plot
- show() - show plot from previous simulation
- disp() - display parameters and initial values from the last simulation
- describe() - describe culture, broth, parameters, variables with values/units


Note that both disp() and describe() takes values from the last simulation and the command process_diagram() brings up the main configuration

Brief information about a command by help(), eg help(simu)
Key system information is listed with the command system_info()

```
plt.rcParams['figure.figsize'] = [30/2.54, 24/2.54]
```

```
# The process diagram is made outside Modelica to illustrate the configuration
process_diagram()
```

No processDiagram.png file in the FMU, but try the file on disk.

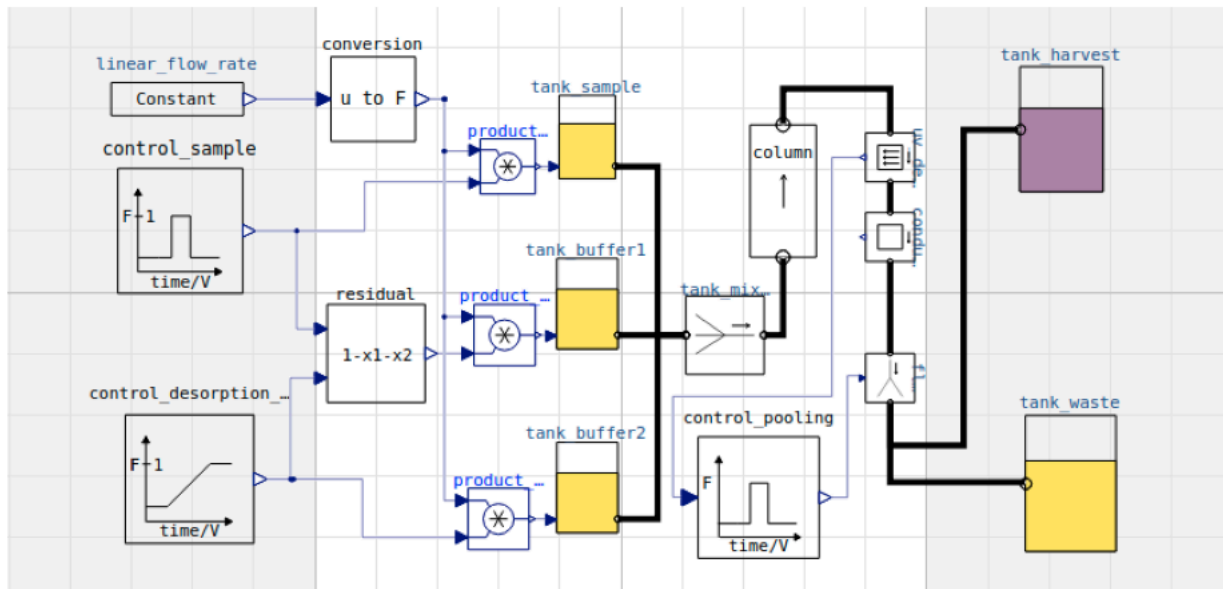

```
describe('chromatography'); # print(); describe('liquidphase')
```

Ion exchange chromatography controlled with varying salt-concentration. The pH is kept constant.

✓ Loading or adsorption

The parameter notation and values are the same as in the referred report. However the flow rate is here denoted F while q in the report. The column is divided in n=8 sections and set at compilation time. The values are arbitrarily chosen in the report and the focus is on qualitative aspects of the model.

The simplified model describe only the column in terms of volume and does not distinguish a high column with a small diameter from a lower with larger diameter.

The parameters k_1 , k_2 , k_3 , k_4 and Q_{av} are given relative volume and with increased column volume a larger capacity is thus obtained.

Loading of the column - try to reproduce Jonas figure 13.

```
newplot(title='Adsorption along the column during loading', plotType='Loading')
```

Sample

```
par(P_in=1.0, A_in=1.0, E_in=0)
```

Column properties

```
par(k1=0.3, k2=0.05, k3=0.05, k4=0.3, Q_av=3.0)
```

```
par(height=20, diameter=0.714)
```

```
par(x_m=0.3)
```

Operation

```
par(E_in_desorption_buffer=8)
```

```
par(LFR=12)
```

```
#par(scale_volume=False)
```

```
par(start_adsorption=0, stop_adsorption=50)
```


```
par(start_desorption=150, stationary_desorption=450)
```

```
par(start_pooling=220, stop_pooling=450)
```

Simulation

```
simu(100)
```

```
↳ Could not find cannot import name 'dopri5' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo/lib/__i
Could not find cannot import name 'rodas' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo/lib/__i
Could not find cannot import name 'odassl' from 'assimulo.lib' (/usr/local/lib/python3.11/site-packages/assimulo/lib/__i
Could not find ODEPACK functions.
Could not find RADAR5
Could not find GLIMDA.
```

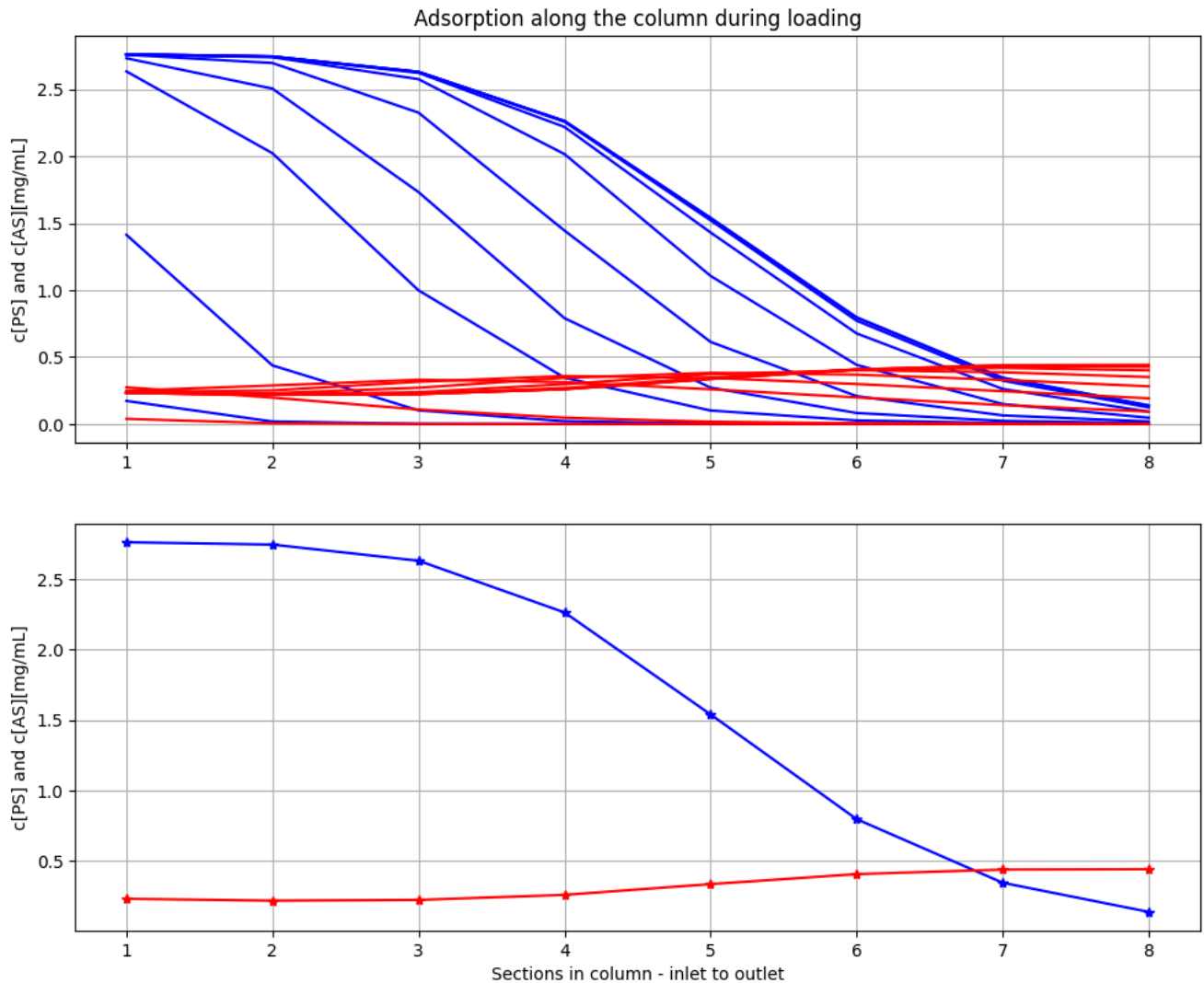

The results are the same as Figure 13 in [1].

We just check that we had the same volume flow rate as Jonas

```
describe('F')
```

```
↳ : 0.08
```

```
describe('V')
```

```
↗ Column volume total - derived : 8.008 [ mL ]
```

```
model.get('column.x_m')
```

```
↗ array([0.3])
```

```
model.get('column.V_m')
```

```
↗ array([2.40235705])
```

```
#describe('column.n')
```


```
model.get('column.column_section[1].V_m')
```

```
↗ array([0.30029463])
```

```
# Impact of shorter time for loading and then less material
newplot(title='Loading', plotType='Loading-combined')
show()
```

```
# Simulation with changed parameter t2
par(stop_adsorption=30); simu(100)
```

```
# Reset changed parameter
par(stop_adsorption=50)
```

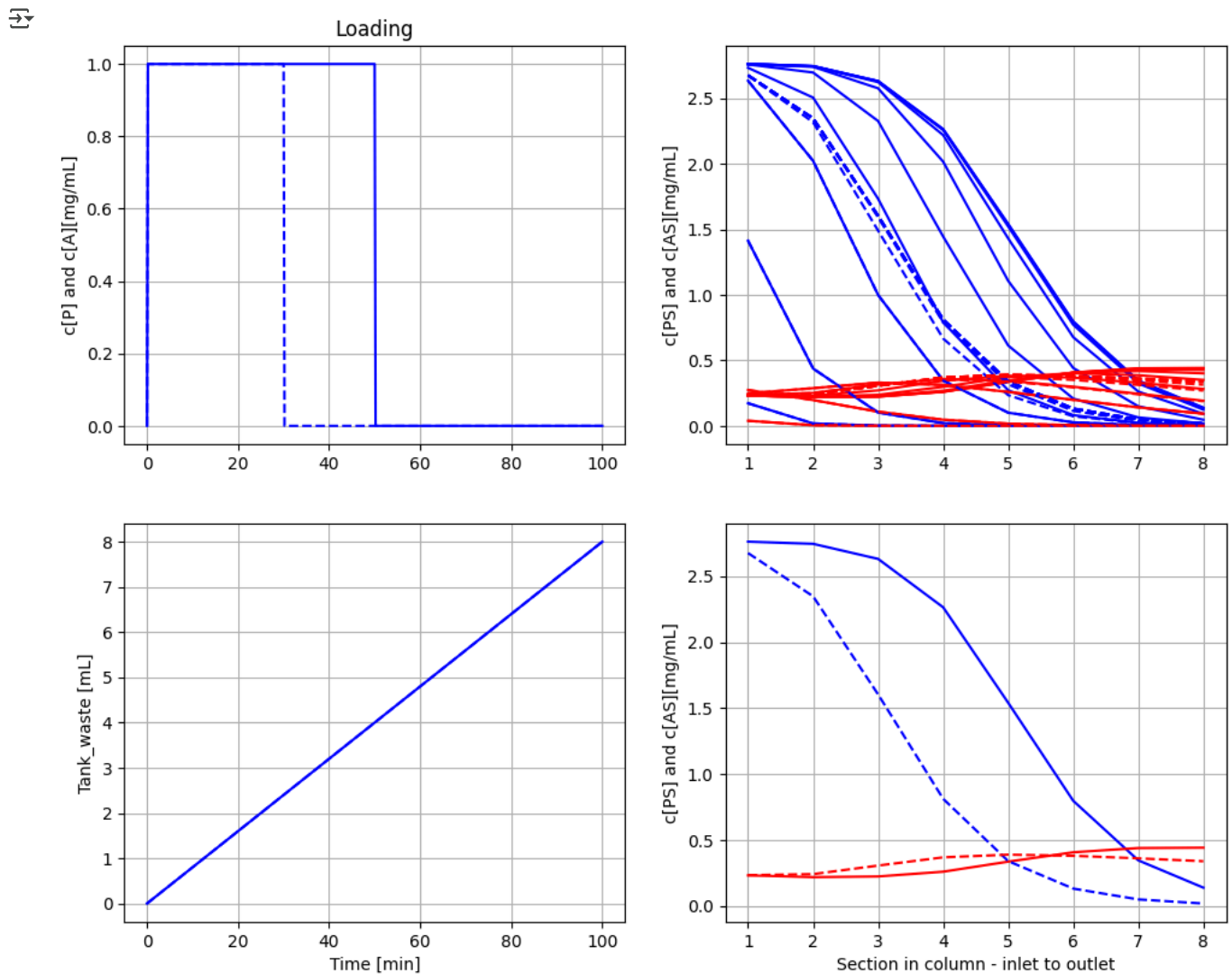

To the left the inlet loading over time. To the right upper concentration along the column at different times and in steady states finally To the right lower concentrations along the column in steady state.

We see that a shorter time and then less material makes less of the column capacity used.

Note that the flow through the column is constant despite change from sample to just buffer 1, and shown in how the volume of the waste tank increase with time.


✓ Elution or desorption

```
# Elution of the column
newplot(title='Elution peaks', plotType='Elution')

# Sample
par(P_in=1, A_in=1.0, E_in=0)

# Operation
par(E_in_desorption_buffer=8)
par(LFR=12.0, start_adsorption=0, stop_adsorption=50, start_desorption=150, stationary_desorption=450)

# Simulation
simu(550)
```

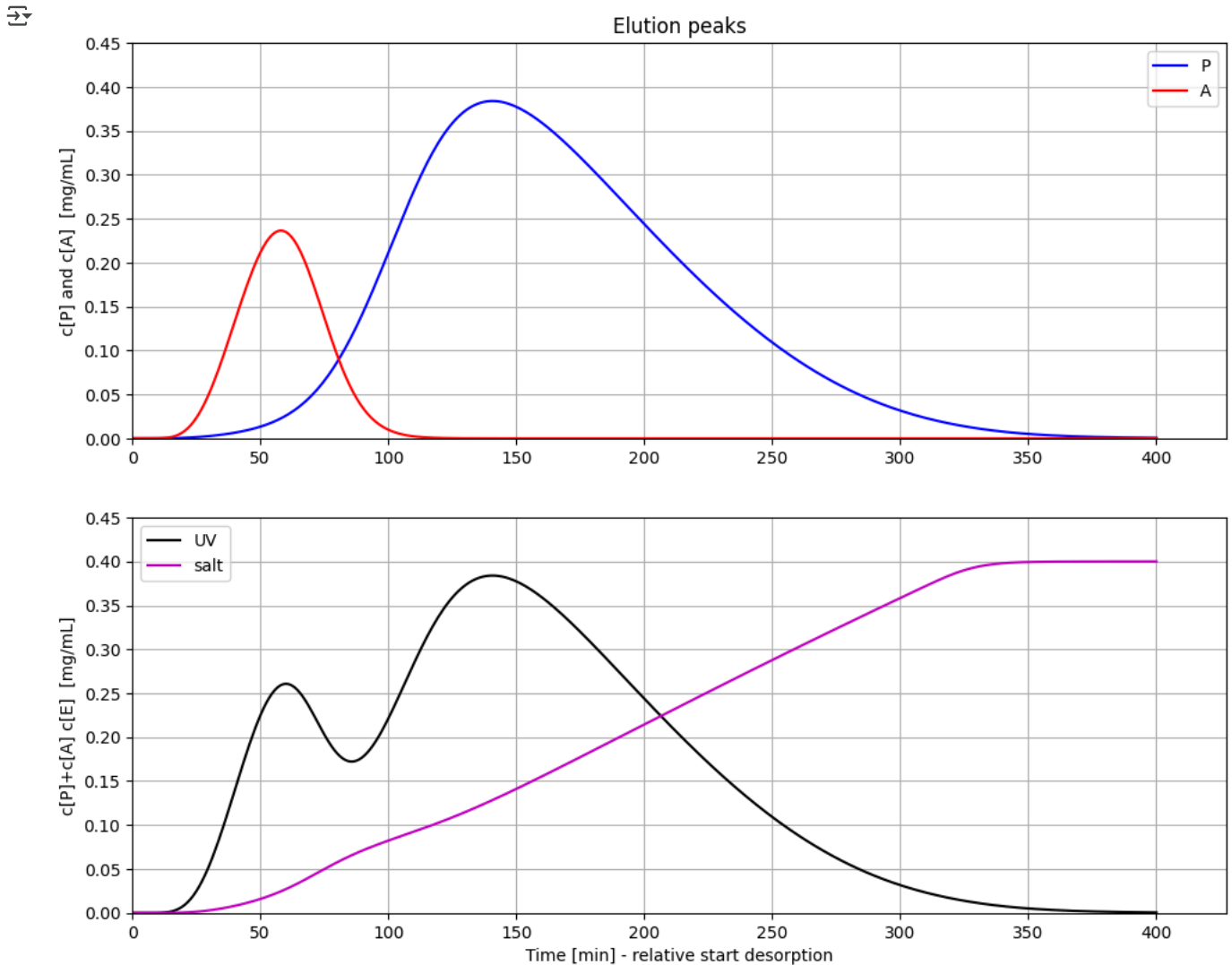

The results are the same as Figure 14 in [1].

The upper diagram shows the column outlet concentrations of P and A over time.

The lower diagram shows the sum (or possibly the UV signal) at column outlet as well as the salt concentration. We have some separation between the two peaks.

Note that the salt concentration deviates slightly from the linear increase between time 50 to 100. This is due to ion interaction with P and A in the column. This phenomenon can also be seen in real data. The ion-salt concentration is scaled with factor 0.05 to get comparable concentrations to P and A.

I have here simulated time 150 of adsorption and then started elution. Here time is counted as zero at the time of start of elutions. Not sure how long Jonas simulated to get steady state before he did elution.

```
# More complete visualization of the elution phase and the different flows
newplot(title='Elution and corresponding flows for two different loadings', plotType='Elution-combined')
```


```

par(stop_adsorption=50); simu(550)

# Simulation with changed parameter t2
par(stop_adsorption=30); simu(550)

# Reset changed parameter
par(stop_adsorption=50)

```

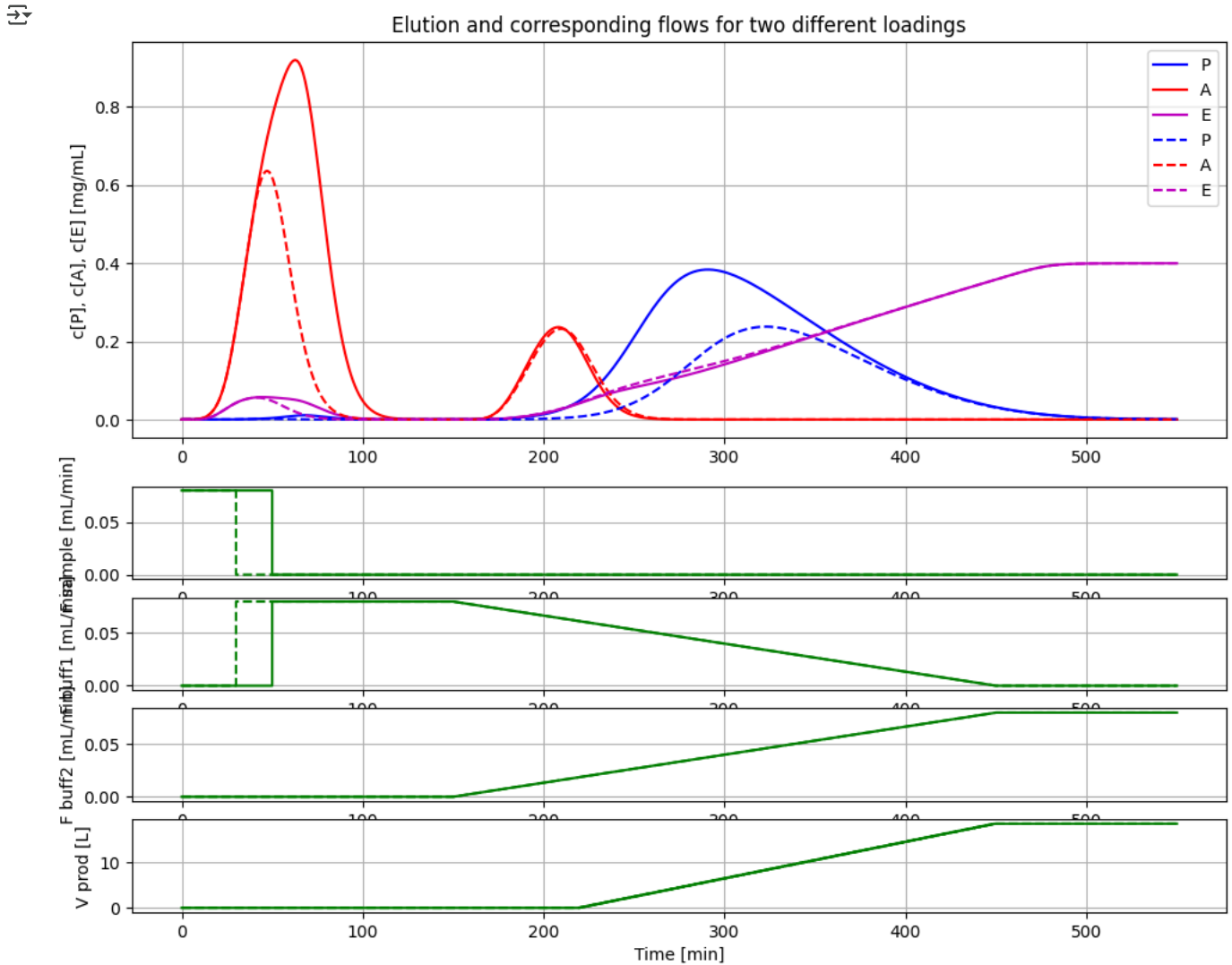

Here a diagram that shows the peaks at the outlet as shown in the previous diagram. Below the flow rates of the three different sources. Here time is 0 at start of adsorption and elution starts at time 150.

Automatic pooling based on UV-measurement is tested in another notebook.

✓ Change of resin properties

```

# Elution of the column
newplot(title='Elution peaks for different resin parameters k1 and k2', plotType='Elution')

# Sample
par(P_in=1, A_in=0.0, E_in=0)


# Operation
par(E_in_desorption_buffer=8)
par(LFR=12.0, start_adsorption=0, stop_adsorption=50, start_desorption=150, stationary_desorption=450)

# Simulations
par(k1=0.05, k2=0.50); simu(550)
par(k1=0.05, k2=0.25); simu(550)
par(k1=0.05, k2=0.05); simu(550)
par(k1=0.25, k2=0.05); simu(550)
par(k1=0.50, k2=0.05); simu(550)

# Adjust diagrams
ax1.set_ylim([0, 0.65])

```

```
ax2.set_ylim([0, 0.65])
plt.show()
```

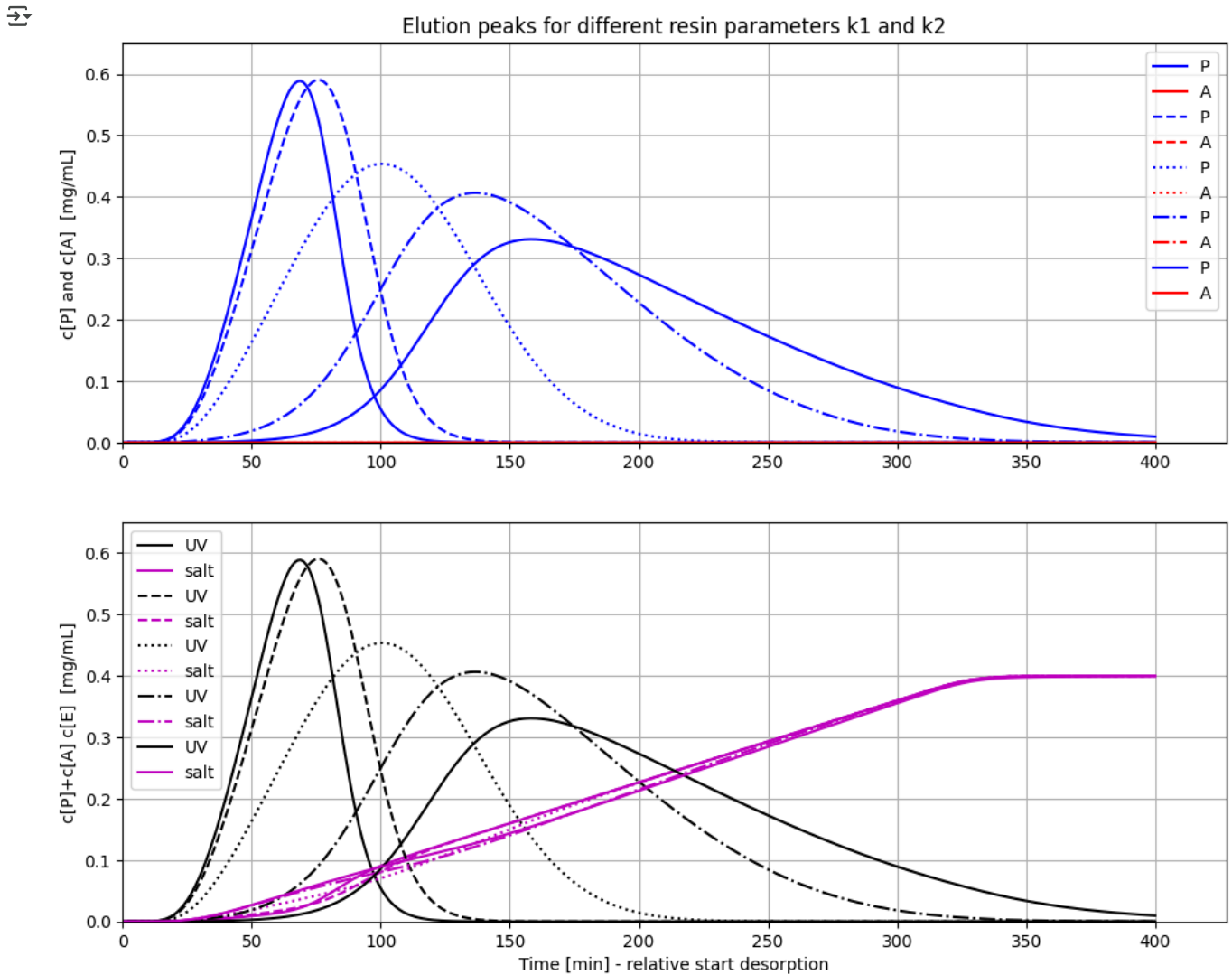

The results are the same as Figure 17 in [1].

Summary

Three important diagrams Figure 13, 14 and 17 in the original report [1] were reproduced and the implementation in Modelica used here is considered validated.

The model is now extended with an improved parametrization of the column that matches industrial practice.

Acknowledgement


The author thanks Karl Johan Brink for sharing his know-how of chromatography operation. He has especially given input to how to parametrize the model in terms often used in industrial practice.

References

- 1) Månsson, Jonas, "Control of chromatography column in production scale", Master thesis TFRT-5599, Department of Automatic Control, LTH, Lund Sweden, 1998.
- 2) Pharmacia LKB Biotechnology. "Ion Exchange chromatography. Principles and Methods.", 3rd edition, 1991.

✓ Appendix

```
describe('MSL')
```


 MSL: 3.2.3 – used components: RealInput, RealOutput, CombiTimeTable, Types

system_info()



System information
-OS: Linux
-Python: 3.11.11
-Scipy: not installed in the notebook
-PyFMI: 2.16.3
-FMU by: OpenModelica Compiler OpenModelica 1.25.0~dev-133-ga5470be
-FMI: 2.0
-Type: FMUModelME2
-Name: Column_system
-Generated: 2024-11-07T11:26:29Z
-MSL: 3.2.3
-Description: Bioprocess Library version 2.3.0
-Interaction: FMU-explore version 1.0.0

Start coding or [generate](#) with AI.